


Mass Spectrum List Report

Calibrante Formiato de sodio - Modo negativo

Acquisition Parameter

Source Type	ESI	Ion Polarity	Negative	Set Nebulizer	2.0 Bar
Focus	Not active	Set Capillary	3500 V	Set Dry Heater	200 °C
Scan Begin	100 m/z	Set End Plate Offset	-500 V	Set Dry Gas	7.0 l/min
Scan End	1000 m/z	Set Collision Cell RF	250.0 Vpp	Set Divert Valve	Source

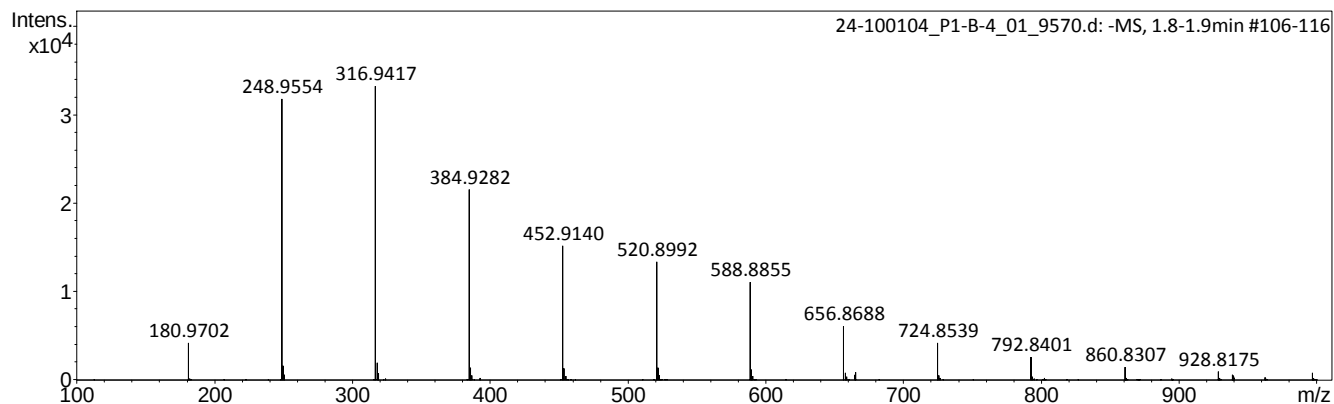

#	m/z	I	I %
1	180.9702	4180	12.6
2	248.9554	31767	95.7
3	249.9581	1553	4.7
4	250.9592	546	1.6
5	316.9417	33186	100.0
6	317.9449	1971	5.9
7	318.9444	741	2.2
8	384.9282	21471	64.7
9	385.9300	1368	4.1
10	386.9322	539	1.6
11	452.9140	15146	45.6
12	453.9179	1285	3.9
13	520.8992	13316	40.1
14	521.9024	1382	4.2
15	588.8855	11015	33.2
16	589.8880	1141	3.4
17	656.8688	6031	18.2
18	657.8711	783	2.4
19	664.8650	525	1.6
20	664.9633	556	1.7
21	665.5760	836	2.5
22	665.6621	732	2.2
23	724.8539	4188	12.6
24	725.8604	537	1.6
25	792.8401	2609	7.9
26	860.8307	1455	4.4
27	928.8175	967	2.9
28	939.1861	599	1.8
29	939.3524	504	1.5
30	996.7993	821	2.5