

Mass Spectrum List Report

Analysis Info

Analysis Name D:\Data\ggc\EF4-S_P1-D-3_01_9676.d
 Method lc_ms directos tunelow300924.m
 Sample Name EF4-S
 Comment Sv: H2O + MeOH + HCOOH
 D'Accorso


Acquisition Date 11/8/2024 3:54:14 PM

Operator EB - MV

Instrument / Ser# microTOF-Q

Acquisition Parameter

Source Type	ESI	Ion Polarity	Positive	Set Nebulizer	2.0 Bar
Focus	Not active	Set Capillary	4000 V	Set Dry Heater	200 °C
Scan Begin	100 m/z	Set End Plate Offset	-500 V	Set Dry Gas	7.0 l/min
Scan End	1000 m/z	Set Collision Cell RF	250.0 Vpp	Set Divert Valve	Source

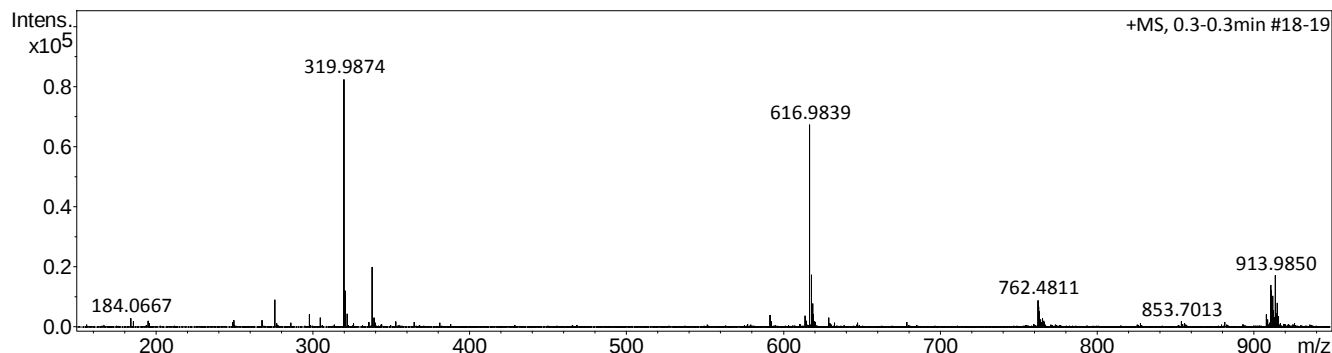

#	m/z	I	I %
1	184.0667	2868	3.5
2	276.0243	9042	11.0
3	298.0048	4142	5.0
4	305.0787	3033	3.7
5	319.9874	82111	100.0
6	320.9901	12066	14.7
7	321.9852	4316	5.3
8	337.9978	19894	24.2
9	339.0009	3058	3.7
10	591.4801	3918	4.8
11	613.9785	3731	4.5
12	616.9839	67079	81.7
13	617.9861	17362	21.1
14	618.9839	7711	9.4
15	619.5044	4124	5.0
16	628.9747	3096	3.8
17	762.4811	8882	10.8
18	762.9800	6614	8.1
19	763.4776	5292	6.4
20	765.4854	2852	3.5
21	907.9634	4076	5.0
22	910.9764	13975	17.0
23	911.4808	11905	14.5
24	911.9791	10287	12.5
25	912.4751	5423	6.6
26	912.9778	3336	4.1
27	913.9850	17081	20.8
28	914.4835	4291	5.2
29	914.9886	7879	9.6
30	915.9829	3485	4.2