

Qualitative Compound Report

Analysis Info


Analysis Name D:\Data\ggc\24-022703_P1-E-2_01_8653.d
 Method lc_qtof_directospos.m
 Sample Name 24-022703
 Comment FAM248
 Sv: MeOH + HCOOH
 Florencia Musso

Acquisition Date 3/4/2024 4:21:57 PM

Operator ELB
 Instrument micrOTOF-Q II

Acquisition Parameter

Source Type	ESI	Ion Polarity	Positive	Set Nebulizer	3.5 Bar
Focus	Not active	Set Capillary	4000 V	Set Dry Heater	200 °C
Scan Begin	70 m/z	Set End Plate Offset	-500 V	Set Dry Gas	7.0 l/min
Scan End	1000 m/z	Set Collision Cell RF	300.0 Vpp	Set Divert Valve	Source

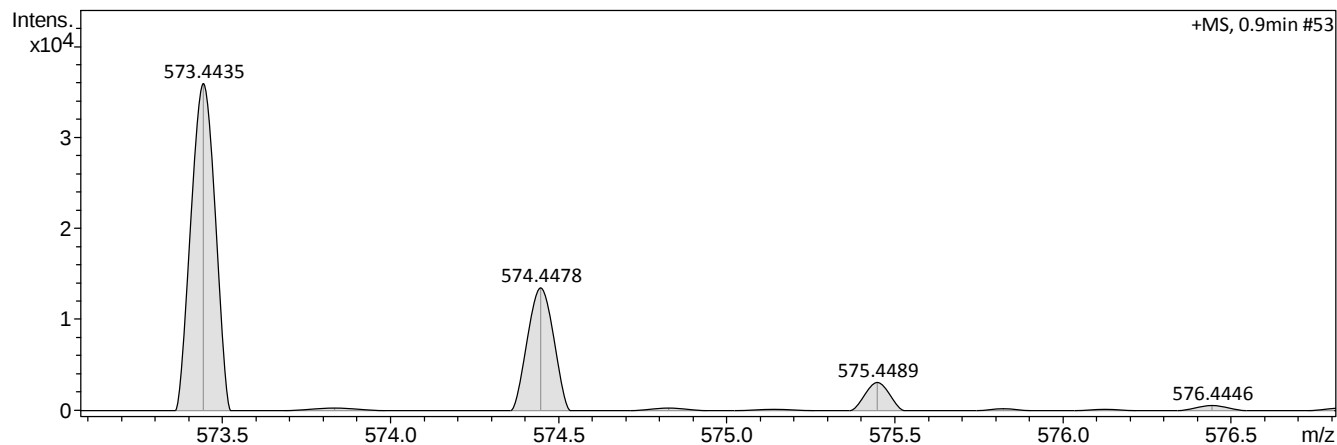

Meas. m/z	#	Ion Formula	m/z	err [ppm]	mSigma	# mSigma	Score	rdb	e ⁻ Conf	N-Rule
573.4435	1	C38H57N2O2	573.4415	-3.6	26.1	2	50.66	11.5	even	ok

Qualitative Compound Report

Analysis Info


Analysis Name D:\Data\ggc\24-022704_P1-E-3_01_8654.d
 Method lc_qtof_directospos.m
 Sample Name 24-022704
 Comment FAM249
 Sv: MeOH + HCOOH
 Florencia Musso

Acquisition Date 3/4/2024 4:27:03 PM

Operator ELB
 Instrument micrOTOF-Q II

Acquisition Parameter

Source Type	ESI	Ion Polarity	Positive	Set Nebulizer	3.5 Bar
Focus	Not active	Set Capillary	4000 V	Set Dry Heater	200 °C
Scan Begin	70 m/z	Set End Plate Offset	-500 V	Set Dry Gas	7.0 l/min
Scan End	1000 m/z	Set Collision Cell RF	300.0 Vpp	Set Divert Valve	Source

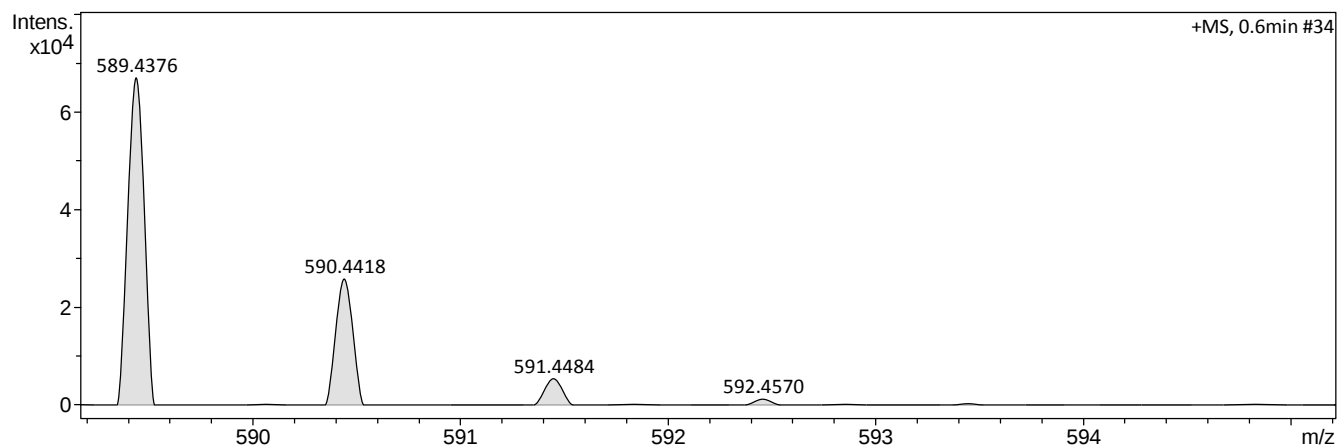

Meas. m/z	#	Ion Formula	m/z	err [ppm]	mSigma	# mSigma	Score	rdB	e ⁻ Conf	N-Rule
589.4376	1	C38H57N2O3	589.4364	-2.0	21.4	1	100.00	11.5	even	ok

Qualitative Compound Report

Analysis Info


Analysis Name C:\Data\24-022705_P1-E-4_01_8655.d
 Method lc_qtof_directospos.m
 Sample Name 24-022705
 Comment FAM344
 Sv: MeOH + HCOOH
 Florencia Musso

Acquisition Date 4/3/2024 16:32:11

Operator ELB
 Instrument micrOTOF-Q II

Acquisition Parameter

Source Type	ESI	Ion Polarity	Positive	Set Nebulizer	3.5 Bar
Focus	Not active	Set Capillary	4000 V	Set Dry Heater	200 °C
Scan Begin	70 m/z	Set End Plate Offset	-500 V	Set Dry Gas	7.0 l/min
Scan End	1000 m/z	Set Collision Cell RF	300.0 Vpp	Set Divert Valve	Source

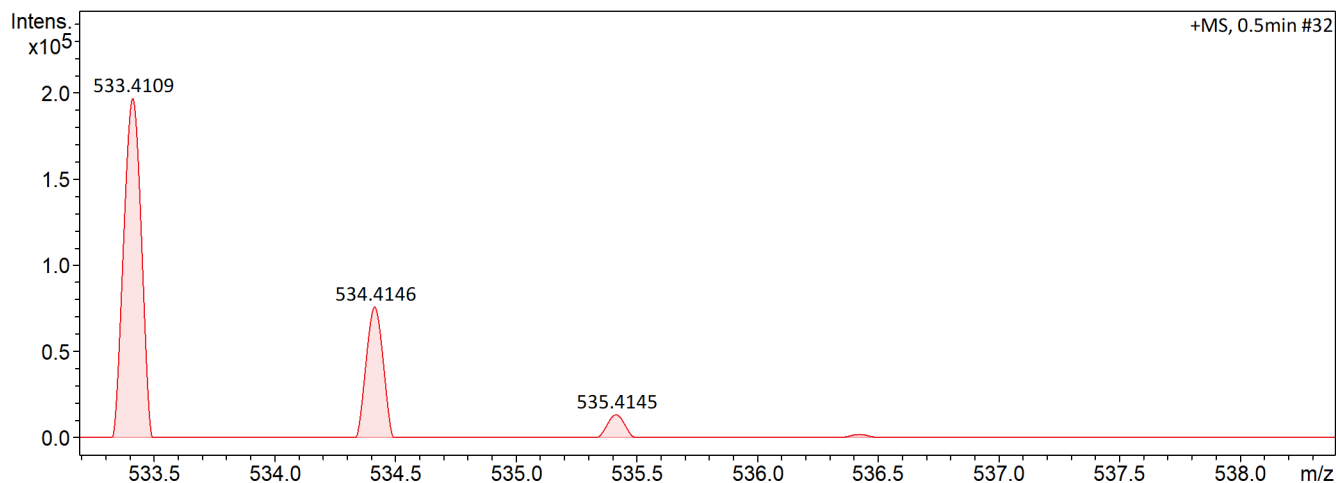

Meas. m/z	#	Ion Formula	m/z	err [ppm]	mSigma	# mSigma	Score	rdb	e ⁻ Conf	N-Rule
533.4109	1	C35H53N2O2	533.4102	-1.3	6.4	1	100.00	10.5	even	ok